

OPERATOR'S MANUAL

.....

APPENDIX C

PACKING PROCEDURE

Introduction

The following is the factory recommended packing procedure for the Infinity laser system. This procedure should be followed if the Infinity system is to be shipped to another location after initial installation.

The Infinity laser system requires two shipping crates, one primarily for the laser head (P/N: 0166-949-00) and the second primarily for the power supply (P/N: 0166-950-00). Table C-1 gives a complete listing of the contents of each crate, in addition to the laser head and power supply, when the system was originally shipped from Coherent.

Table C-1. Contents of Infinity Laser Head and Power Supply Shipping Crates.

Laser Head Crate
1. PC and Laser Documentation
2. Spare Lamp
3. Computer Monitor
4. Deionization Cartridge
5. Computer Keyboard
6. Maintenance Kit
7. Water Hoses

Power Supply Crate
1. Computer (CPU)

Laser Head

The Infinity head crate consists of two compartments, see Figure C-1. The compartment for the Infinity head is foam filled (left side in Figure C-1), it is therefore very important to line the compartment with anti-static material before placing the laser head in the crate, see Figure C-1 and the left hand side of Figure C-2.

To prevent damage to laser head optical components the two POT assembly interlock pins (one on each side of the POT assembly covering the POT clips) should be removed and placed in the maintenance kit. The POT cover over the lamp leads should be secured using tape. Also, tape should be placed over the DPSS output port, and on both ends of the Faraday Isolator. Lens tissue should be taped over the two output ports in the Infinity head cover.

To prevent the umbilical connectors from becoming damaged the end should be enclosed in packing material (e.g., bubble wrap), see the right hand side of Figure C-2.



- 1. Infinity Laser Head
- 2. Anti-Static Wrap

- 3. Second Compartment for Infinity Components

Figure C-1. Overview of Infinity Head Crate



Figure C-2. Infinity Laser Head Properly Packed for Shipment

All other components of the Infinity laser system, with the exception of the computer, should be placed in the second compartment of the Infinity head crate (see Figure C-1, C-3, and Table C-1).



Figure C-3. Packing of Infinity Components, see Table C-1 for Part Descriptions

Power Supply

The Infinity power supply crate also consists of two compartments, see Figure C-4. Like the Infinity head crate, the compartment for the Infinity power supply is foam lined. It is therefore very important to wrap anti-static material around the power supply before placing it in the shipping crate. The front cover of the Power Supply crate should be utilized as a ramp for getting the power supply into the crate, see Figure C-4. The Infinity computer should be placed in the upper compartment of the Power Supply shipping crate. The key for the power supply should be placed in the Infinity Maintenance kit.

Note that the packing foam on the inside face of the Power Supply crate front is removable, see Figure C-5. Removing the piece furthest from the crate, when the cover is installed as a ramp, will allow the Power Supply to be inserted into the crate much more easily.

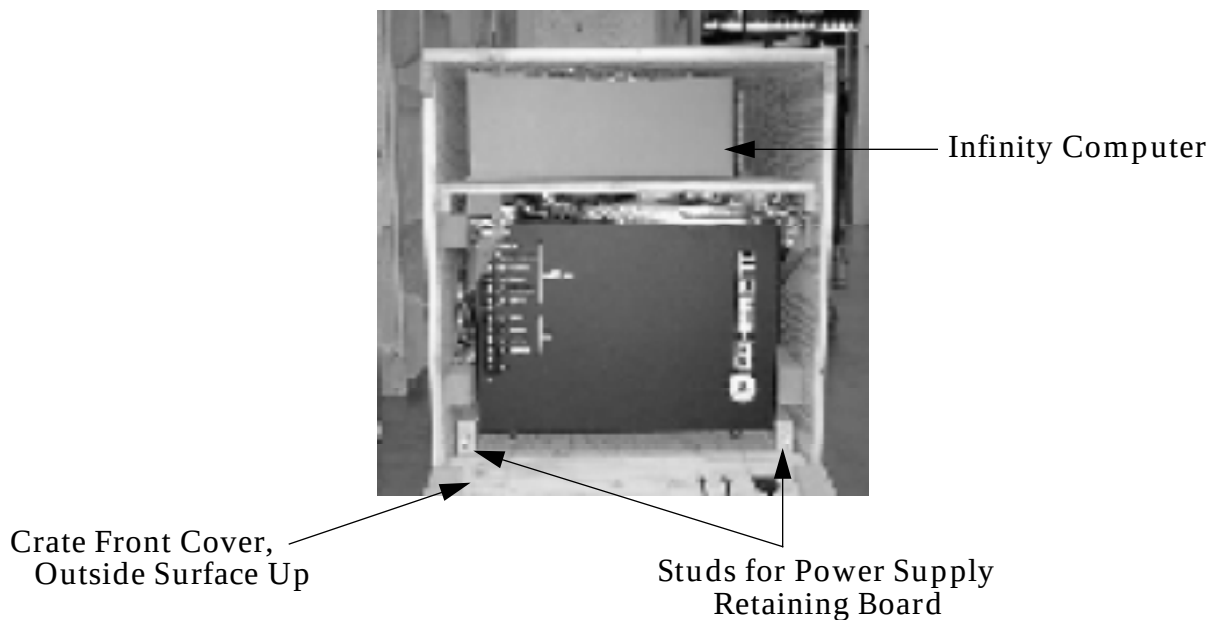


Figure C-4. Infinity Laser Head Properly Packed for Shipment

Once positioned in the crate the power supply is secured using the retaining board and wing nuts as illustrated in Figure C-5.

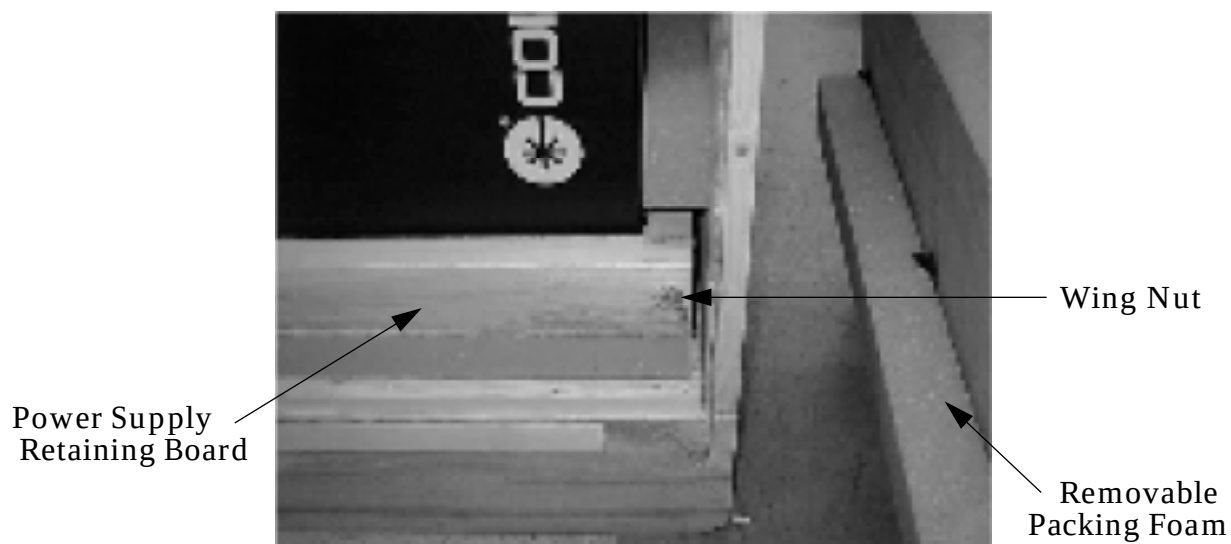


Figure C-5. Installed Power Supply Retaining Board